

How to look after your prism

Cleaning your prism

The rough side of the prism can collect dust and may need to be cleaned.

- Use a gentle stream of cool/tepid water and a very soft brush (e.g. a toothbrush) – clean in the direction of the lines on the prism.
- Pat the prism and glasses dry with a soft cloth
- If the prism becomes dirty with make-up or grease, it may need to be removed and cleaned in the same way with a weak mixture of washing up liquid and water.




Re-fitting your prism

If you need to remove the prism or it falls off, it can be replaced easily

- Wash the glasses lens and the prism with a weak mixture of washing-up liquid or soap to remove any grease. Rinse well with water.
- Submerge the glasses and the prism in a bowl of cool/tepid water or under a stream of running water (1).
- Slide the prism with the smooth side against the back of the glasses lens (2).
- Remove the glasses from the water, holding the prism in place.
- Gently squeeze any water or air bubbles out from between the prism and the glasses lens with your thumb (3).
- The prism will slide easily around the lens at this stage if changes of it's position are needed.
- Pat the prism and your glasses dry using a clean, dry cloth or paper towel.
- It may take a while for the prism to dry onto the lens: be careful with your glasses during this time.



(1)

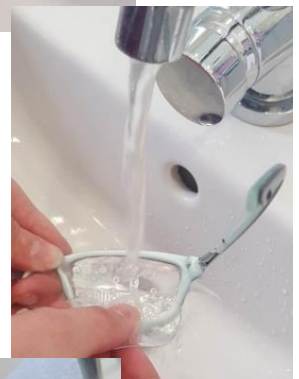

(2)



(3)